MANN+HUMMEL

PRODUCT TECHNICAL BULLETIN



SINGLE-LAYER

Disposable Face Masks

- Lightweight and breathable
- White, single-layer, spunbound material — PP
- Comfortable, elastic strap keeps mask in place
- Compact, foldable design — just pop open and it's ready to use



Color: White

Size: One size fits most adults

Material: Single-layer, Spunbound, Polypropylene (PP), 107 gsm

Head Band: 1/4″ elastic band


MADE IN THE USA

MANN+HUMMEL



Recommendations for use:

- Disposable material — not to be washed for reuse
- Replace mask every 4 to 6 hours for optimum performance
- Not to be worn during sleep
- Store in a well-ventilated room, free of corrosive gas
- Not recommended for use:
  - In any surgical setting or where exposure to liquid, bodily or other hazardous fluids may be expected
  - In clinical settings, where infection risk level through inhalation exposure is high
  - In the presence of high-intensity heat source or flammable gas
  - In places of poor ventilation or where breathing is impeded
  - With children under the age of 3

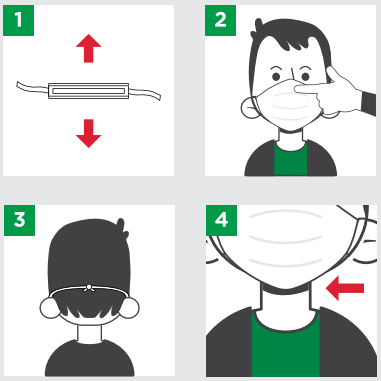

MADE IN THE USA

WASH YOUR HANDS BEFORE HANDLING THE MASK.

1. Hold the mask on the top and bottom and pop it open quickly with a little force.
2. Once open, pull both sides to expand. Position the mask over the face, covering the nose and mouth.
3. Pull elastic up and over your head.
4. Pull down the edge of the mask, so that it fully covers the chin.

FDA Emergency Use Authorization dated 4/20/2018 Required Notice:

- The product has not been FDA cleared or approved.
- The product has been authorized by FDA under an EUA for use by HCP as PPE to help prevent the spread of infection or illness in healthcare settings and by the general public to help slow the spread of the virus during the COVID19 pandemic.
- This product is authorized only for the duration of the declaration that circumstances exist justifying the authorization of the emergency use of medical devices, including alternative products used as medical devices, during the COVID-19 outbreak, under section 564(b)(1) of the Act, 21 U.S.C. § 360bbb-3(b)(1) unless the authorization is terminated or revoked sooner.

©2020 MANN+HUMMEL Purolator Filters LLC
Fayetteville, NC 28306, USA
All rights reserved.
Techline 800-526-4250